



BMP15 ELISA Kit

(For Human)

Enzyme Immunoassay for the quantification of human BMP15 in serum, plasma and cell culture supernatant

Catalog number: ARG80971

For research use only. Not for use in diagnostic procedures.

TABLE OF CONTENTS

SECTION	Page
INTRODUCTION	3
PRINCIPLE OF THE ASSAY	3
MATERIALS PROVIDED & STORAGE INFORMATION	4
MATERIALS REQUIRED BUT NOT PROVIDED	4
TECHNICAL HINTS AND PRECAUTIONS	5
SAMPLE COLLECTION & STORAGE INFORMATION	5
REAGENT PREPARATION.....	6
ASSAY PROCEDURE	6
CALCULATION OF RESULTS	7
EXAMPLE OF TYPICAL STANDARD CURVE	8
QUALITY ASSURANCE.....	8

MANUFACTURED BY:

Arigo Biolaboratories Corporation

Address: No. 22, Ln. 227, Gongyuan Rd., Hsinchu City 300, Taiwan

Phone: +886 (3) 562 1738

Fax: +886 (3) 561 3008

Email: info@arigobio.com

INTRODUCTION

BMP15 protein is a member of the bone morphogenetic protein family which is part of the transforming growth factor-beta superfamily. The transforming growth factor-beta superfamily includes large families of growth and differentiation factors. It is thought that this protein may be involved in oocyte maturation and follicular development as a homodimer or by forming heterodimers with a related protein, Gdf9. Defects in this gene are the cause of ovarian dysgenesis 2. [provided by RefSeq, Sep 2009]

PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for human BMP15 has been pre-coated onto a microtiter plate. Standards or samples are pipetted into the wells and any BMP15 present is bound by the immobilized antibody. A HRP-conjugated antibody specific for BMP15 is added to each well and incubate. After washing away any unbound antibody-enzyme reagent, a chromogenic substrate solution is added to the wells and color develops in proportion to the amount of BMP15 bound in the initial step. The color development is stopped by the addition of acid and the intensity of the color is measured at a wavelength of 450 nm \pm 2 nm. The concentration of BMP15 in the sample is then determined by comparing the O.D of samples to the standard curve.

BMP15 ELISA Kit ARG80971

MATERIALS PROVIDED & STORAGE INFORMATION

Store the unopened kit at 2-8 °C. Use the kit before expiration date.

Component	Quantity	Storage information
Antibody-coated microplate	96 wells	4°C
Standards A-F (Concentration: 62.5, 125,250,500,1000,2000 pg/ml)	6 vials X 0.5ml	4°C
Sample Diluent	6 ml	4°C
HRP-conjugated Antibody	10 ml	4°C
20X Wash buffer	25 ml	4°C
Chromogen solution A	6 ml	4°C (Protect from light)
Chromogen solution B	6 ml	4°C (Protect from light)
Plate sealer	2	RT
STOP solution	6 ml	4°C

MATERIALS REQUIRED BUT NOT PROVIDED

- Microplate reader capable of measuring absorbance at 450nm
- 37 °C incubator
- Pipettes and pipette tips
- Deionized or distilled water
- Automated microplate washer (optional)

TECHNICAL HINTS AND PRECAUTIONS

- Wear protective gloves, clothing, eye, and face protection especially while handling blood or body fluid samples.
- Store the kit at 4°C at all times.
- Briefly spin down the antibody conjugate concentrate and HRP-Streptavidin concentrate before use.
- If crystals are observed in the 20X Wash buffer, warm to RT (not more than 50°C) until the crystals are completely dissolved.
- Ensure complete reconstitution and dilution of reagents prior to use.
- It is highly recommended that the standards, samples and controls be assayed in duplicates.
- Change pipette tips between the addition of different reagent or samples.

SAMPLE COLLECTION & STORAGE INFORMATION

The sample collection and storage conditions listed below are intended as general guidelines. Sample stability has not been evaluated.

Cell Culture Supernatants - Remove particulates by centrifugation and aliquot & store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

Serum - Use a serum separator tube (SST) and allow samples to clot for 30 minutes before centrifugation for 15 minutes at 1000 x g. Remove serum and assay immediately or aliquot and store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

Plasma - Collect plasma using EDTA or heparin as an anticoagulant. Centrifuge for 15 minutes at 1000 x g within 30 minutes of collection. Assay immediately or aliquot and store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

REAGENT PREPARATION

- **1X Wash buffer:** Dilute 20X Wash buffer into distilled water to yield 1X Wash buffer.
- **Samples:** If samples give OD readings higher than the highest standard provided, please dilute samples with sample diluent and repeat the assay.

ASSAY PROCEDURE

All materials should be equilibrated to room temperature (RT) before use. Standards and samples should be assayed in duplicates.

1. Remove excess microplate strips from the plate frame, return them to the foil pouch containing the desiccant pack, and reseal it.
2. Add 50 μ l Sample Diluent Buffer in duplicate into wells (blank).
3. Add 50 μ l of standards and samples in duplicate into wells.
4. Add 100ul of HRP-conjugated antibody to each well.
5. Mix well gently, cover and incubate the plate for 1 hour at 37°C.
6. Aspirate each well and wash, repeating the process 4 times for a total 5 washes. Wash by filling each well with 1 \times Wash Buffer (350 μ l) using a squirt bottle, manifold dispenser, or autowasher. Complete removal of liquid at each is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating, decanting or blotting against clean paper towels.
7. Add 50 μ l of Chromogen solution A and 50 μ l of Chromogen solution B to each well subsequently. Cover wells and incubate for 15 minutes at 37°C. Protect from light.

BMP15 ELISA Kit ARG80971


8. Add 50 μ l of Stop Solution to each well. Mix well.
9. Read the OD with a microplate reader at 450 nm immediately.

CALCULATION OF RESULTS

1. Calculate the average absorbance values for each set of standards and samples.
2. Using linear graph paper, construct a standard curve by plotting the mean absorbance obtained from each standard against its concentration with absorbance value on the vertical (Y) axis and concentration on the horizontal (X) axis.
3. Using the mean absorbance value for each sample determine the corresponding concentration from the standard curve.
4. Automated method: The results in the IFU have been calculated automatically using a 4 PL (4 Parameter Logistics) curve fit. 4 Parameter Logistics is the preferred method. Other data reduction functions may give slightly different results.

EXAMPLE OF TYPICAL STANDARD CURVE

The following data is for demonstration only and cannot be used in place of data generations at the time of assay.

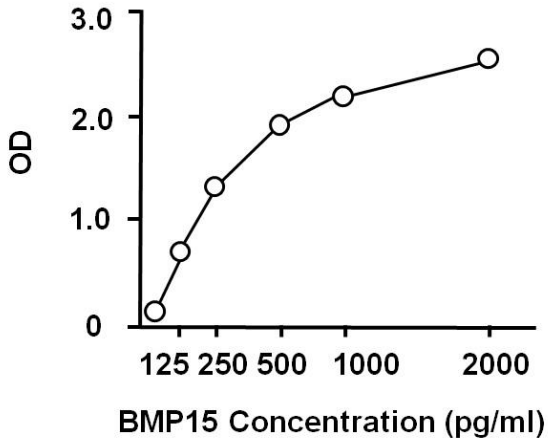

QUALITY ASSURANCE

Sensitivity

The minimum detectable dose (MDD) of BMP15 ranged from 62.5-2000 pg/ml.

The mean MDD was 10 pg/ml.